USER MANUAL

Massage gun 23575

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Pistolet do masażu 23575

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Massagepistole 23575

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Pistolet de massage 23575

INSTRUCCIONES DE OPERACION

pistola de masaje 23575

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Pistola da massaggio 23575

MANUAL DE UTILIZARE

Pistol de masaj 23575

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Masážna pištoľ 23575

KASUTUSJUHEND

Massaažipüstol 23575

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Masažo pistoletas 23575

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Masāžas pistole 23575

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
00023575 (M31)

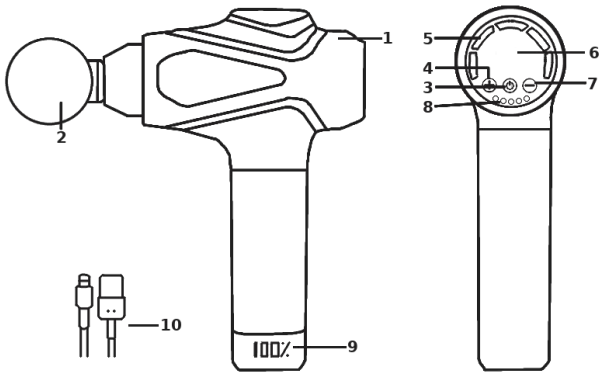EN

**Ladies and Gentlemen, thank you for purchasing our product!
Before starting to use, please read the following instructions
for the correct use of the product. Please keep these instructions
on
future and to follow its recommendations, as non-compliance may endanger life or health.**

Important: KEEP FOR FUTURE REFERENCE!

- Please read it completely before use and follow all instructions in the manual.
- Do not use the product if any part is missing or damaged.
- Please keep this manual for future reference.

SPECIFICATION:

- Black colour.
- Power: 45W.
- Motor rated voltage: 12V.
- Battery type: lithium battery (Li-Ion).
- Capacity: 2500 mah.
- Charging time: 120 min.
- Working time: 6 h.
- Dimensions (length/width/height): 25 x 18 x 7cm.
- Cable length: 100 cm.

CHARACTERISTIC:

The massager allows continuous operation, up to 6 hours on one charge. Its performance is at the highest level. Take care of proper muscle regeneration - the gun will allow you to massage at home, without the need for constant visits to a physiotherapist. It will relieve muscle soreness, stiffness and pain. It has as many as 30 levels of massage intensity. This is a very large range that allows you to individually adjust the speed to the muscles in which it is to be used. The gun has a clear display that shows a given level, so you have constant control. The intensity of work is changed using a sensitive touch panel, just touch it gently - you don't have to force it to change the mode. Each muscle group requires appropriate stimuli, and thanks to 10 different tips, it is possible to work with many muscles, even those that are difficult to access. The device is made of the best quality materials. It fits perfectly in the hand, the place where we hold it is covered with an anti-slip structure that prevents it from slipping during the procedure. The set is packed in a practical case/suitcase with a handle, and inside there is an insert with special embossing for each of the set's elements. On the other side, there is a mesh compartment that will allow you to pack e.g. training materials. The case will allow you to take the set with you when traveling, and it will be easy to transport for mobile physiotherapists.

INCLUDED:

- Case.
- Pistol.
- USB cable.
- Charger.
- 12 tips.

PART DESCRIPTION (SEE PAGE 2):

- 1) Massage gun body
- 2) Massage head
- 3) Power button
- 4) Gear up button
- 5) Battery usage percentage
- 6) Display screen
- 7) Down gear key
- 8) Pressure sensor
- 9) Battery indicator
- 10) Charging cable
- 11) Ball massage head - suitable for muscle groups such as buttocks and thighs.
- 12) U-shaped massage head - suitable for shoulder spine massage.
- 13) D-shaped massage head - For relaxing and shaping muscles in all parts of the body.
- 14) Curved massage head - relaxing massage of the arm and leg muscles.
- 15) Shock-absorbing massage head - suitable for affecting deep tissues such as the meridian joints of the hand.
- 16) Flat heat for massage - relaxation and shaping of muscles in all parts of the body.
- 17) Teeth Massage Head - Apply to neck, abdomen, arms and legs to eliminate fat massaging acupuncture points.
- 18) Tower massage head - suitable for women with small calf, abdominal and arm muscles.
- 19) Conical massage head - suitable for precise massage of pain points.
- 20) T-shaped massage head - suitable for modeling massage of all legs.
- 21) Cooling head - a soft, gel tip that supports muscle regeneration after strength or running training and the fight against pain of various origins. To prepare the cooling tip for work, place it in the refrigerator for 1-2 hours so that its temperature drops and you can start the cooling massage. It can also be used without prior cooling to massage delicate parts of the body.
- 22) Pneumatic (air) head – used for a very gentle massage of the most tense muscle groups. It is ideal for massaging the bone area and can also be an alternative to the round tip for more sensitive people. It is also great for a relaxing massage.

USER MANUAL:

1) Charging

- Fully charge the battery (up to six hours) before first use: plug the adapter into a wall socket, then connect the USB cable to the adapter and the other end to the battery charging port.
- Once the battery is connected, the LEDs stop flashing and remain fully illuminated.
- Full charge is indicated when all LEDs flash and remain fully illuminated.
- The battery can be charged at any battery level.
- It is not recommended to completely discharge the battery to the level of the red LED.

2) Before use, prepare the skin or thin clothing on the part of the muscle to be relaxed. The massaged surface should be dry and clean.

3) Starting: turn the switch at the bottom of the handle from OFF to ON. To turn the massager on or off, press the power button.

4) Turn off: flip the switch at the bottom of the handle.

HOW TO USE:

1) Install the massage head: According to different parts of use and effects, choose the appropriate massage head. When installing the massage head, please insert the massage head directly into the outlet of the massage gun and insert it firmly all the way so that it is flush with the line of the massage head.

2) Use the massage gun: turn on the gun, switch to the appropriate gear and adjust the appropriate pressure. Gently press the massage head against the muscle group to be relaxed, move it slowly back and forth, then stay there for 30 seconds after finding the sensitive and sore spot. The entire massage process should not exceed 10 minutes. The gun has enough pressure. Avoid bones, joints, arteries and other parts.

MAINTENANCE AND CLEANING

- Pull out the massage head directly, wipe the head and gun body, clean the housing with a slightly damp soft towel, and then use paper to dry it.
- The product should be stored in a dry, well-ventilated and non-corrosive place.
- Do not use corrosive solutions to wipe the case, otherwise it may cause product failure or product cracking, discoloration.
- When transporting, the product should be fixed and placed gently to prevent violent vibration, collision, tipping and dropping.
- When storing or moving, please clean it and put it back into the packing box or storage box.

WARNINGS:

- To reduce the risk of injury, electric shock, fire, equipment damage, or malfunction, read all important information and warnings before use.
- It is strictly prohibited to use this product by injured, sick or minors. In case of muscle injury, consult a doctor before using the device.
- Do not use on wet or dirty skin. The gun can be used in places covered by clothing. The massager should be pressed against the body and gently moved between different points, pausing at each point for about 60 seconds.
- Do not use the massager on parts of the body above the neck or in areas with less muscle or bone. The device may only be used on soft tissues if the user does not experience pain or other adverse reactions.
- Do not use the device on body parts with bruises or if there is pain or unusual body reactions.
- Do not immerse the product in water or allow water to enter the product's ventilation holes.
- Do not keep it near fire or heat sources.
- Do not throw the device and avoid accidental falls.
- Do not disassemble or modify the product.
- Fingers, hair or other body parts may become trapped in the shaft and back of the massage head.
- Children under 8 years of age must not use the device. Children must not use, clean or maintain the appliance without adult supervision.
- People who meet any of the following conditions should consult a doctor before using the massager: pregnant women, people with a pacemaker, people who have recently had surgery, people with spine, joint or bone defects, people who have complications caused by other diseases (such as neuropathy or retinal damage); modern scientific research focuses on determining the impact of vibration massage on people with special health conditions and may limit or expand the scope of the above-mentioned contraindications. Practical experience has shown that the combination of treatment and vibration massage can also serve as a type of therapy in many circumstances. However, such an approach should be determined and ordered by a specialist doctor or therapist or carried out under special supervision in their presence.
- It is forbidden to run this product without the massage head installed.
- Before use, please check whether the product can operate normally according to the operation method of this manual.
- Before each use, carefully check whether the device or tip is undamaged and can be used normally according to the methods specified in the user manual.
- The room temperature of the use environment is 0-40°C.
- The room temperature of the charging environment is 0-40°C.
- The time of use of the same location should not exceed 60 seconds.
- Please clean the device after use to avoid the build-up of residue such as sweat or grease.
- Always unplug the charger after charging is complete.
- Do not use the product while charging.

INFORMATION REGARDING BATTERY AND ACCUMULATION

- Keep batteries/rechargeable batteries out of the reach of children.
- Do not store batteries/rechargeable batteries in easily accessible places, as there is a risk that they may be swallowed by children or pets. If swallowed, consult a doctor as soon as possible.
- Leaking or damaged batteries/rechargeable batteries may cause burns when in contact with the skin, so always wear appropriate protective gloves when touching them.
- Batteries/rechargeable batteries must not be short-circuited, disassembled or thrown into fire. There is a risk of explosion!
- Ordinary disposable batteries should not be charged, as there is a risk of explosion! Only charge batteries intended for this purpose, using an appropriate charger.
- Charge the batteries before inserting them into the product
- Batteries should only be charged under the supervision of an adult.
- Do not mix alkaline, standard (carbon-zinc) and rechargeable (nickel-cadmium) batteries.
- Never mix new and old batteries.
- It is forbidden to interfere with the power supply cables - it may cause a short circuit.
- Exhausted batteries should be removed from the product.
- Do not throw used batteries into unsorted waste!
- Batteries should be disposed of in accordance with the regulations for the disposal of hazardous electronic waste.
- If batteries are not disposed of correctly, hazardous substances may pose a risk to human health or the environment. Recycling materials contributes to the conservation of natural resources.
- Always insert batteries/rechargeable batteries with the polarity correct (plus/+ and minus/-).

Information for consumers on how to handle used electrical and electronic equipment:

Waste electrical and electronic equipment must not be disposed of with other household waste. Such equipment should be collected selectively, as indicated by the attached marking - crossed out wheeled waste container. Failure to comply with the above rule and improper disposal of used equipment may pose a threat to the natural environment and human health resulting from the presence of hazardous components in the equipment, such as: electrical wiring, plastic, switches, etc. To avoid the above threat, these components should be removed, collected and appropriately processed by specialized companies. Disassembling used equipment yourself is unacceptable. Households play an important role in contributing to reuse and recovery, including recycling of waste equipment. This is done in particular through active participation in the waste equipment collection system.

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